4-Stretcher Ambulance





SHORT DESCRIPTION

The WAS 4-stretcher ambulance has been designed for transporting casualties in the event of major catastrophes. The vehicle can transport up to four people, under medical supervision, to a nearby treatment centre for further care. All four patients can therefore be given emergency medical treatment. The stretcher brackets to the left and right can be folded up against the side walls to save space and create extra room. If necessary, the stretcher brackets are lowered in a user-friendly way to enable ergonomic loading.



Spineboards, scoop stretchers or stretcher chairs are stored on the rear doors.



The oxygen cylinders are stowed behind the sliding side door making them easily accessible, and at the same time saving space.



While patients are being transported, the loading ramps are simply folded away, providing additional stability.

The concept is transferable to other base vehicle types. The equipment is an example and can be customized.

EXTERIOR

Configurable warning systems for optimum visibility in emergency situations, all lighting available in LED or halogen options

Audible warning device in the front lightbar or installed separately in the bumper 360° scene lights to illuminate the treatment and surrounding areas

INTERIOR

Fold-away frame on the left and right side walls, each for carrying two emergency stretchers Seamless interior trim for ease of infection control Hi-grip and fully sealed flooring in assorted colour options. Composite floor structures also available High-performance air-conditioning or full climate control unit in patient compartment

Infusion bottle holders for each patient

Two-seater bench with three-point seat belts on

the partition wall for medical staff

Bulkhead drawer unit for storing

medical equipment

Oxygen supply for all four patients

High level wall lockers on the front left side

Storage for stretcher, carry chair, scoop stretcher and spineboards

Overhead storage compartment above the cab,

accessed from the patient compartment

Centre console between driver's seat and co-driver's seat in the cab

Optional: Base vehicle with air conditioning in cab

Optional: Intercom between cab and patient compartment Optional: Reversing camera


Optional: Camera system in the patient compartment

BASE VEHICLE

Model	Mercedes-Benz Sprinter/
	VW Crafter panel van
Wheel base	3665 mm
Total length	6086 mm
Total width	1933 mm
Total height	2630 mm
Maximum	3880 kg
total weight	



The stretcher brackets to the left and right can be folded up against the sidewalls to maximise space.

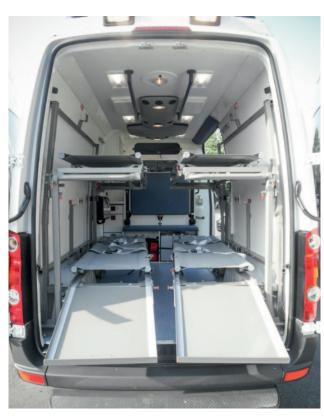

The stretcher brackets can be lowered to enable ergonomic loading.

