

WAS[®]

MAKING VEHICLES SPECIAL



Medical Care On-Site: The **WAS 900 Mobile Medical Practice.**

WAS | 900 MOBILE MEDICAL PRACTICE

The WAS 900 Mobile Medical Practice: The Solution for Mobile Medical Care.



SHORT DESCRIPTION

The changes in demographics and advances in mobile medical treatment are some of the challenges that face the modern medical services. That is the reason why we have made medical treatment mobile. The mobile medical practice is a fully equipped doctors practice on wheels. It is suitable for various medical fields, for the use at large social events (instead of emergency ambulances) and for medical pre- and post-care. With a total vehicle weight of 3.5 t, it can be driven with a normal car driver's license. Due to its modular furniture the mobile medical practice can be easily adjusted to the sensitive demands regarding the diagnosis and medical care of refugees. The separated dressing room enlarges the comfort zone for the patients. The modular box body is crash tested according to EN 1789 and is equipped with modern, energy efficient LED lighting, doctor's spots and blue calming trauma lights. An independent vehicle heater and 230-V-airconditioning ensure pleasant temperatures during all seasons. The box body is highly energy efficient due to excellent insulation.



The dressing room is spacious and can be closed with a curtain for privacy.



The moveable surgery bed enables easy access to the patient from all sides. Energy efficient LED lighting ensures perfect conditions.



The lighting can be switched to relaxing blue light.

INTERIOR

Modular interior for different medical applications

Surgery bed with different fixation points for flexible patient access

Rapid exchangeable drawer units for fast adaption to different requirements

Refrigerator for medical supplies

Dressing room


Separate treatment area

Motor independent vehicle heater

230-V-air conditioning

Wash basin with hot and cold water

20 l fresh and waste water tanks

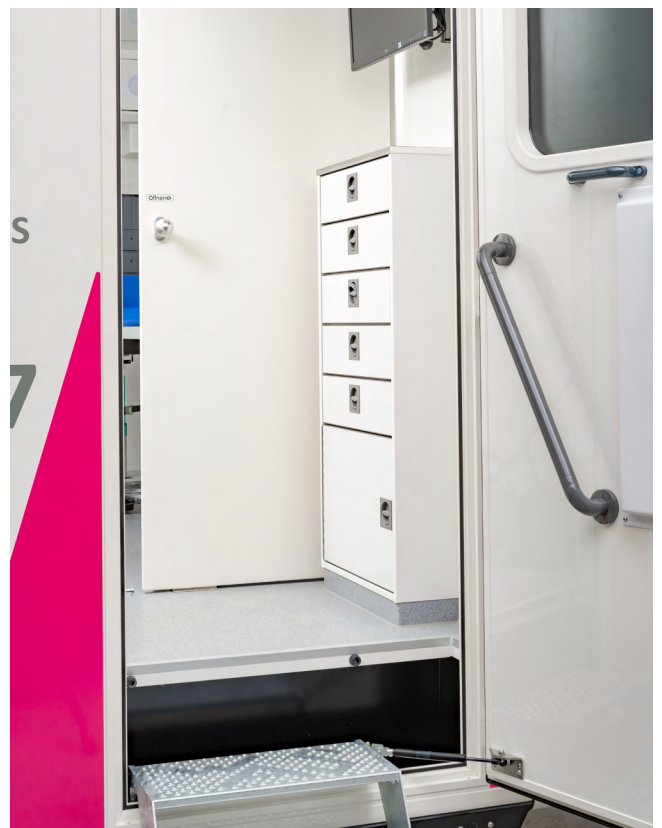

Access to the patient compartment is via a separable entrance area at the rear of the vehicle.

The Mobile Medical Practice: Highly Flexible Due to Modular Interior Design.

DRIVER'S COMPARTMENT

Individual seats in the driver's cab, 4way adjustable, incl. seat heater

Front, side, head- and shoulder air bags

Lane assist and fatigue warner

Radio-CD-player with IVECO SYNC incl. emergency assistant

SPECIAL FEATURES

Crash tested box body construction according to EN 1789 (ambulance standard)

Highly energy efficient due to excellent insulation properties

3.5 t gross vehicle weight, can be driven on a car licence

Modern, energy efficient LED lighting

Roof spoiler

Roof suitable for antenna installation

Optional: Fixed installation of a satellite antenna (Starlink) so that digital services can be accessed anywhere

BASE VEHICLE

Model	Iveco Daily
Type	3,0 I DEVe
Engine power	129 kW (176 PS)
Gearbox	H I - MATIC 8-speed automatic transmission
Wheel base	3450 mm
Maximum total weight	3500 kg



The wash basin can be easily covered to make additional work surface.

The concept is transferable to other base vehicle types. The equipment is an example and can be customized.

